

sunbather



UNCHAINED
LABS

Jump to GMP in a jiff

Taking LNPs and other nanomedicines into GMP can be a costly, painful process with re-validation of different equipment, flow paths and mixers. Sunbather makes the move to GMP a walk in the park with buttoned-up documentation, replaceable flow path and easy transfer built right in. It uses the same mixing process and protocol as its Sunny Suite family members Sunscreen and Sunshine, so you can take your dialed in formulation into GMP without all the fuss.


- Up to 1.8 L/h
- Integrated in-line dilution
- Fully replaceable flow path
- Range of mixer types
- 21 CFR Part 11 compliant software
- Transferrable process from Sunshine



Mix it up big time

Prep one input channel with payloads, the other with lipids or other nanoparticle building blocks, and you're ready to mix up your nanoparticles with a Sunny. GMP-ready Sunnies have identical mixing designs to those used on Sunscreen and Sunshine, so your nailed down process from development sticks the landing into GMP.

Sunny Trident

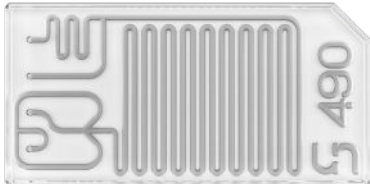


- Reverse-angle mixing
- In-line dilution

Sunny X



- Cross-type mixing
- Multiple channel sizes

Sunny T

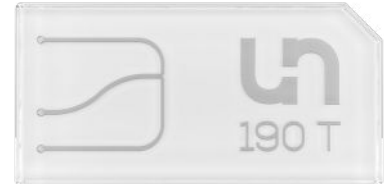

- T-mixing
- Multiple channel sizes

GMP-ready when you are

Sunbather comes available with an installation qualification and operational qualification (IQ/OQ), and comprehensive user training so that you can get up and running as fast as possible. With fully traceable wetted materials, extractables data, and a manufacturing process grounded in by our Quality Management System (QMS), you'll be ready to dive straight into GMP with Sunbather.



No-nonsense, compliant software

Sunbather’s software serves up Administrator and Operator access levels, a 21 CFR Part 11 compliant audit trail, and password signoff for process critical actions. The software’s simple user interface means fast familiarization, and the touchscreen kiosk keeps things clean and makes setting up a new run a breeze.

UNCHAINED
LABS

Sunbather

Run Parameters

	Enable Channel	<input checked="" type="radio"/> mL/min <input type="radio"/> µL/min	Input Volume	Current State
Channel 1	<input checked="" type="checkbox"/>	<input type="text" value="7.5"/> mL/min	187.5 mL	Idle
Channel 2	<input checked="" type="checkbox"/>	<input type="text" value="2.5"/> mL/min	62.5 mL	Idle
Channel 3	<input checked="" type="checkbox"/>	<input type="text" value="10"/> mL/min	250 mL	Idle
Channel 4	<input type="checkbox"/>	<input type="text" value="Channel 4 flow rate"/>	0 mL	Not Initialized
Total Volume		<input type="text" value="500"/> mL		
Total Time		<input type="text" value="0:25:00"/>		
Batch ID		<input type="text" value="UL_Test_Batch"/>		

CFR

Log Out

Power Off

Initialize

Empty

Fill

Stop

START

Device Status :

About

The Sunny Suite

The Sunny Suite has a killer solution for every step of your nanoparticle development.

Sunscreen automates high-throughput screening of up to 96 formulations, then **Sunshine** homes in on the perfect mixing process conditions, with both automated experiments and continuous flow modes. When you're ready for clinical trials, Sunbather dials up your particles in a fully GMP-compliant environment.



Screen




Optimize



Produce

Trust the process

With identical fluid delivery and microfluidic mixing at every stage, the Sunny Suite makes for smooth sailing in every transition from formulation screening to process development to GMP production. Whether you're screening through hundreds of formulations or producing liters of particles for clinical trials, the Sunny Suite's the total package.

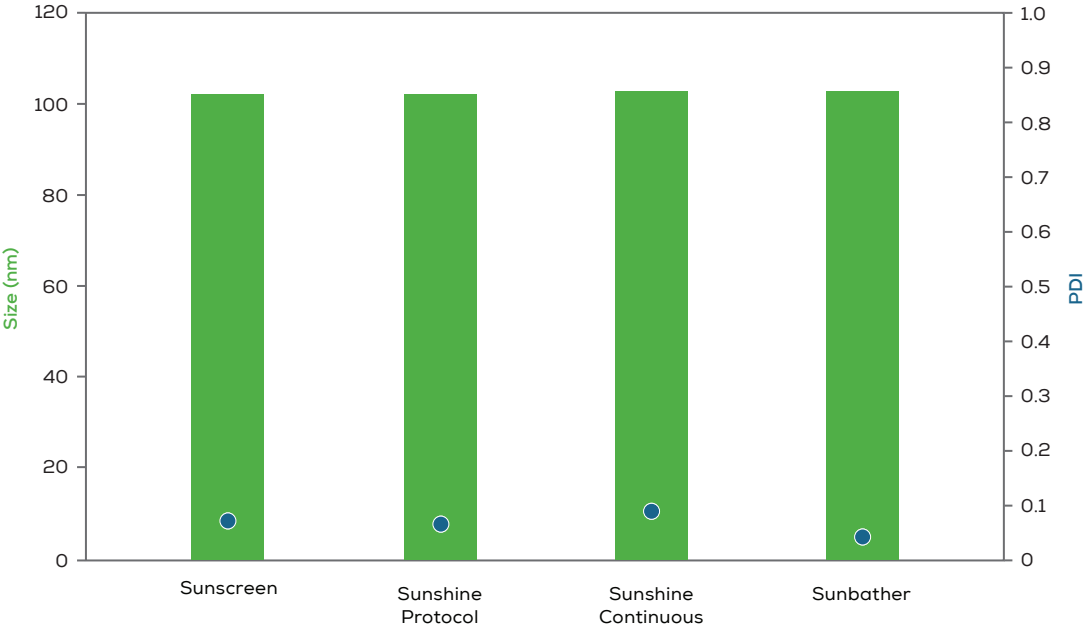

Specifications

Application	
Throughput	10 mL/min per channel
Total flow rate range	0.1–30 mL/min (configuration dependent)
Flow rate ratio (aqueous to organic) range	1:1 to 10:1
Typical particle size range	40–200 nm*
PDI	<0.2*
Encapsulation efficiency	>90%*
Instrument	
Physical	
Weight and dimensions (instrument only)	19 kg, 29 cm H x 41 cm W x 32 cm D
Operating pressure	0–10 bar
Electrical	
Voltage Input	100 V–240 V AC, 50–60 Hz
Power consumption (peak/average)	1A/500 mA
Other information	
Fluidic port sizes	¼"–28 threaded nut, luer fitting (male and female)
Wetted materials	PEEK, FEP, Glass, PTFE, Silicone
Computer	Windows 11 Enterprise Touchscreen PC with Kiosk-mode operation
Consumables	
GMP Sunnies	GMP-ready glass microfluidic mixing devices of various geometries and channel sizes
Available types	Sunny 490 Trident T, Sunny 100 X, Sunny 100 T, Sunny 190 X, Sunny 190 T, Sunny 275 X, Sunny 275 T, Sunny 150 3D, Sunny 50 Micromixer
Sunbather 1-Sunny Flowpath Pack	A complete wetted flowpath for use with the Sunbather instrument running a 1-Sunny configuration.
Sunbather 2-Sunny Flowpath Pack	A complete wetted flowpath for use with the Sunbather instrument running a 2-Sunny configuration.

* Formulation dependent



Unchained Labs
4747 Willow Road
Pleasanton, CA 94588
Phone: 1.925.587.9800
Toll-free: 1.800.815.6384
Email: info@unchainedlabs.com

© 2025 Unchained Labs. All rights reserved. The Unchained Labs logo, Sunbather and Sunbather logo are trademarks and/or registered trademarks of Unchained Labs. All other brands or product names mentioned are trademarks owned by their respective organizations.

Rev B